


# PosiSep® X INTRANASAL SPLINT



- Biocompatible-made of all natural Chitosan-based polymers
- Minimizes bleeding and edema after surgery
- Easily conforms to nasal anatomy
- Offers easier placement with a superior field of view
- No special storage conditions required

**PosiSep® X** is indicated for use as an intranasal splint intended to minimize bleeding and edema and to prevent adhesions between the septum and the nasal cavity. **PosiSep® X** is placed in the nasal cavity after surgery or trauma. The **PosiSep®** line is constructed of a patient-comfortable sponge which is manufactured from naturally occurring Chitosan polymers. Subsequent irrigation will cause the product to fragment, thus optimizing patient tolerance and comfort.

- Fragments by way of gentle irrigation and aspiration
- Easy to apply-no mixing or additional prep time necessary
- Widely studied and published anti-adhesion properties of N,O-carboxymethyl chitosan (NOCC)

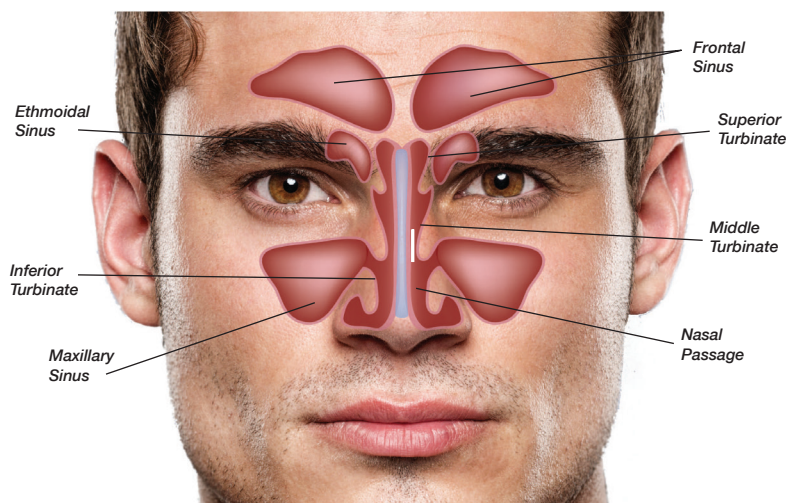

*\*diagram not to scale\**

A small, smooth PosiSep® X sponge is inserted into the nasal cavity where it begins to quickly expand after irrigation. PosiSep® X can be cut into smaller sections for optimal placement.

## Before Hydration



## After Hydration



When PosiSep X is hydrated, it expands to fill the space preventing post-op adhesions.

## What is Chitosan?

Chitosan (ky-toe-san) comes from a substance known as chitin, an amino-polysaccharide removed from the powdered shells of crustaceans.

**CHITOSAN TECHNOLOGY**

**ORDERING INFORMATION** PosiSep® X is available in two sizes. See ordering information below.

### Part Number:

### Description:

9210564

PosiSep® X Intranasal Splint 0.5" x 1.5" 5pk (US)

9210584

PosiSep® X Intranasal Splint 0.6" x 2.0" 5pk (US)

### How Supplied:

PosiSep® X is packaged five (5) sterile, single use units per carton. One (1) sterile unit contains one (1) Intranasal Splint. PosiSep® X is terminally sterilized by gamma irradiation.

CAUTION: Federal (U.S.) law restricts this device to sale by or on the order of a physician. PosiSep® is contraindicated for patients who have known allergies to shellfish. For detailed information regarding indications for use, warnings and precautions, see Instructions for Use.

**To place an order or for more information on PosiSep® or our other medical technologies, please contact a sales representative by calling 866-612-2568.**